

Urine Measurement Systems

Ureofix[®] 500 Classic




Introduction



The **Urefix® 500 Classic** is a well known Urine Measurement System dedicated to the diuresis and the evaluation of the quality of urines within a day. It allows an accurate urine volume measurement. The rigid 3 compartments container has a capacity of 500 ml.

The new generation product benefits from new technologies and fulfils B. Braun environmental standards.


**DEHP FREE
LATEX FREE**

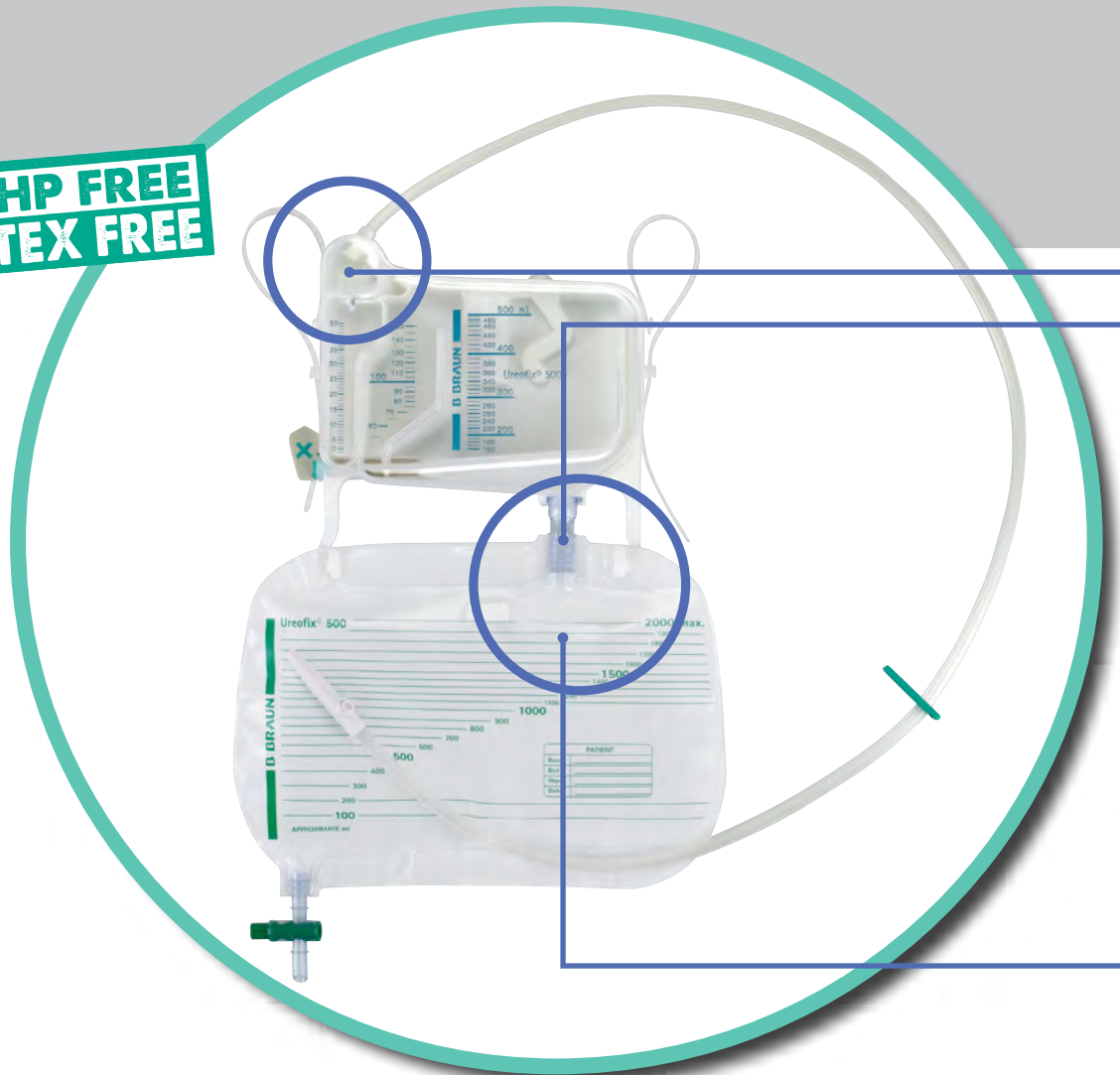

Urefix® 500 Classic - 2 l


Premounted 2 l bag, to facilitate users' quality of life

- **Container drip chamber block**

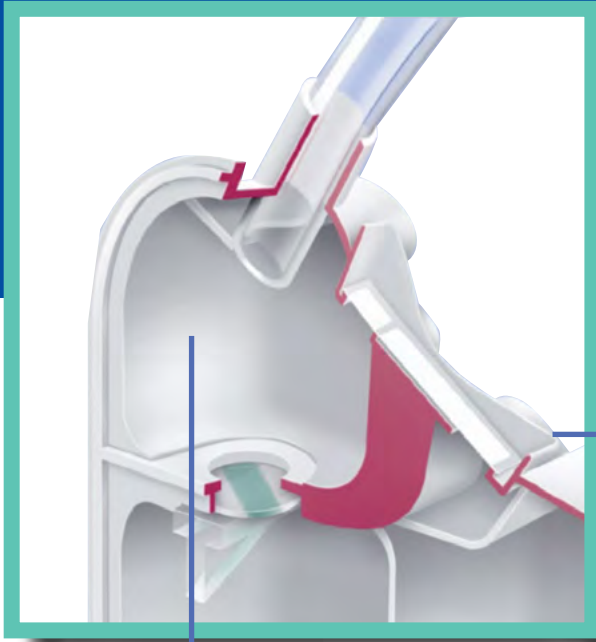


Venting air filters

- Prevent retrograde migration of bacterias
- An integrated Anti Reflux Valve to prevent the retrograde flow of urine into the tubing
- Pressure balance provided thanks to the venting air inlet filters

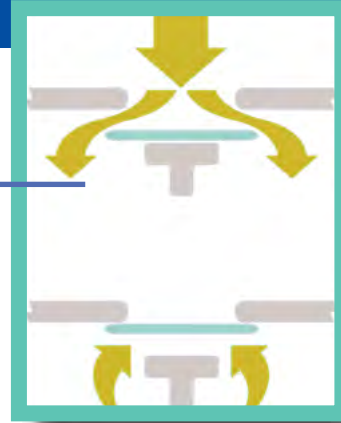

Anti Reflux Valve

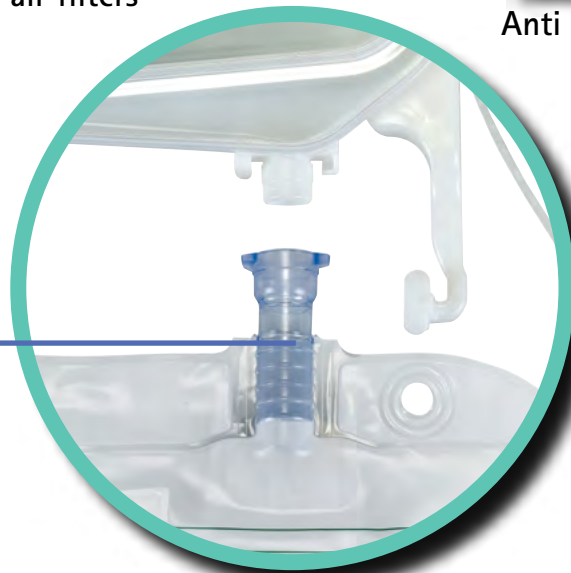

- **Simple twist bayonet connector**
 - An easier and safe connection without leakage
 - Available in 1.5 l / 2 l and 3.5 l



- **Bag with Anti-Reflux Valve and air inlet**

To avoid the clamp usage
Available on 1.5 l, 2 l and 3.5 l bags

- **New bag design with clear patient identification space**



Socket

- A large choice of fixation solutions for the Ureofix® 500 Classic to be adapted on all types of beds and supports

2 integrated fixations

- Integrated removable hook premounted on the back side of the container
- 2 fixation straps premounted on the container



2 additional fixations in option



- Pin and tip jack system



- Removable fixation plate




Mini ifu strap fixation

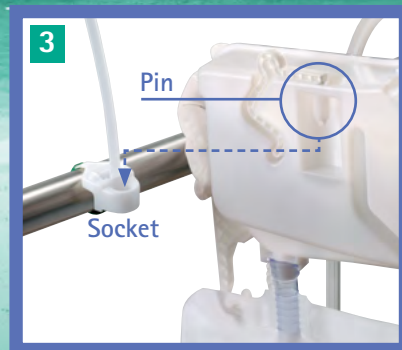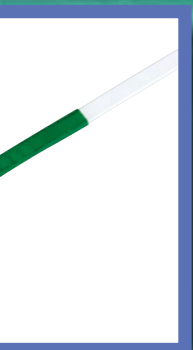

User Benefits

- Reduced risk of infection
- Precise observation
- Easy and safe handling
- Fast emptying of the container
- Safe sampling procedure to avoid needle stick injuries

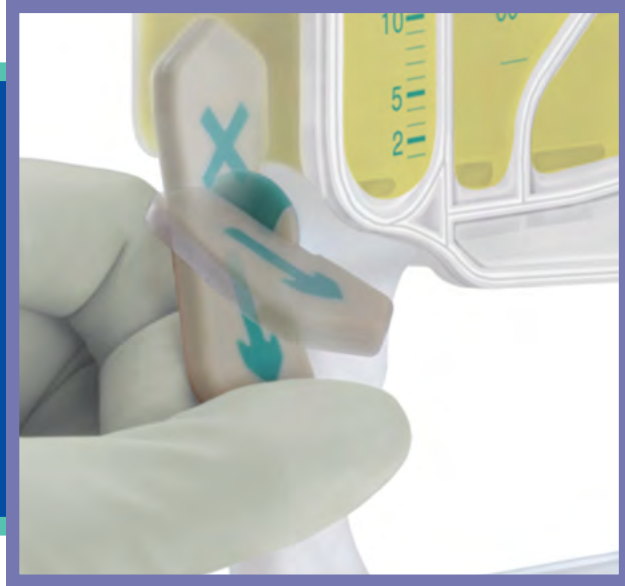

Easy identification of the flow direction

- **A needle free sample port**

The whole range of the new Ureofix 500 Classic includes a needle free sample port function on the tubing connector to match the latest guidelines for a better safety of the user.

- **Sample port compatible with needle free syringes**



- **Horizontal storage of cross-outlet**
to decrease the risk of infection

Product Description and References

Ureofix® 500 Classic Sets

References	Product description	Tubing Lght cm	units/box
With closed bags			
441 7910	Ureofix® 500 Classic with exchange bag 1,5 l	120	10
441 7920	Ureofix® 500 Classic with exchange bag 1,5 l	170	10
441 7960	Ureofix® 500 Classic with exchange bag 3,5 l	170	10
With bottom outlet bags			
441 7930	Ureofix® 500 Classic with premounted b.o bag 2 l	120	10
441 7940	Ureofix® 500 Classic with premounted b.o bag 2 l	170	10
441 7950	Ureofix® 500 Classic with premounted b.o bag 2 l	200	10

Ureofix® 500 Classic Spare Bags

References	Product description	Size	units/box
441 7535	Ureofix® 500 Classic closed bag	1.5 l	100
441 7543	Ureofix® 500 Classic closed bag	3.5 l	100
441 7551	Ureofix® 500 Classic Bottom Outlet Bag	2 l	100

Ureofix® 500 Classic Accessories

References	Product description	units/box
441 7755	Fixation strap	10*10

EO sterilization process

B. Braun Medical S.A.S. | RCS Nanterre 562050856 | 204, avenue du Maréchal Juin
92660 Boulogne Billancourt | France
www.bbraun.com